

VISTA® 637 LSO

Doctor: _____ Fitter: _____

Patient Name: _____ Date: _____

Patient #: _____ Additional Follow-Up Dates: _____

TOOLS NECESSARY: Scissors • Heat Gun • Tape Measure

FOR USE WITH PRODUCTS MANUFACTURED BY ASPEN MEDICAL PRODUCTS ONLY. THIS PRODUCT IS INTENDED FOR APPLICATION BY HEALTH CARE PRACTITIONERS AS DIRECTED BY A PHYSICIAN OR OTHER QUALIFIED MEDICAL AUTHORITY. THIS IS A PREFABRICATED ORTHOSIS. IT IS INTENDED TO BE CUSTOMIZED TO AN INDIVIDUAL PATIENT. FOLLOW THE STEPS BELOW TO CUSTOMIZE.


STEP 1 - MEASUREMENTS

1 Lower Rib Circumference = _____

2 Hip Circumference = _____


3 T9 to Sacrococcygeal Junction = _____

TIME SPENT: _____

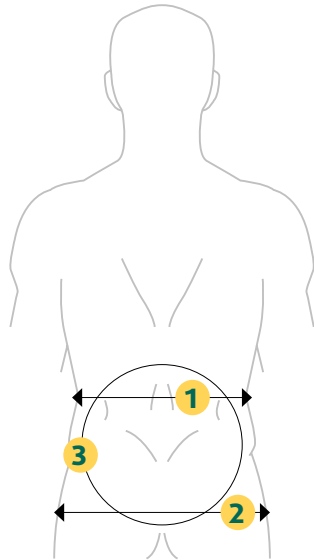

STEP 2 - CUSTOMIZE BACK PANEL TO ANATOMY

- A. Measure patient's lordosis then customize back panel to anatomy.
- B. To customize back panel, remove the panel, heat, trim, and reassemble.

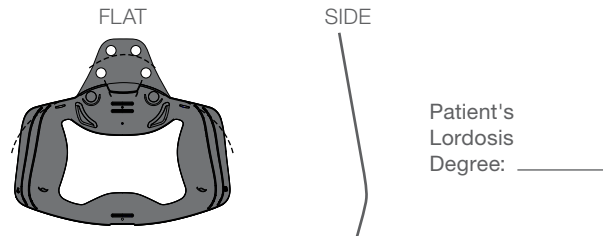

Heat form individual patient's anatomy and contour to create intimate fit for individual's lordosis and soft tissue. Trim for individual patient's anatomy based on 3 _____


- C. Remove lordotic pad to accommodate for lordosis. YES NO

TIME SPENT: _____

STEP 3 - CUSTOMIZE SIZING AND TIGHTENING MECHANISM


SIZING IS CRITICAL TO PROPER PERFORMANCE
Use the measurements below to customize to patient's anatomy.

- A. Use waist circumference (average of 1 and 2 _____) to determine where to fit rivets of belt through proper sizing adjustment holes in sides of back panel.

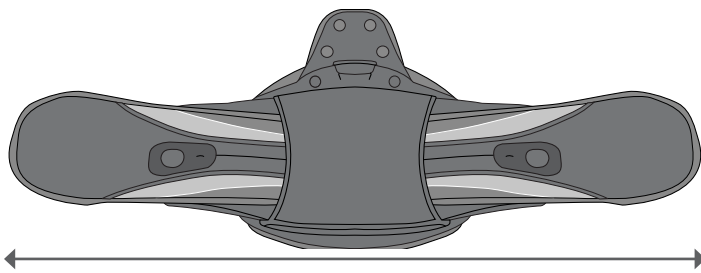

A. _____


- B. Once proper size is achieved, pull taut to lock rivets in place.

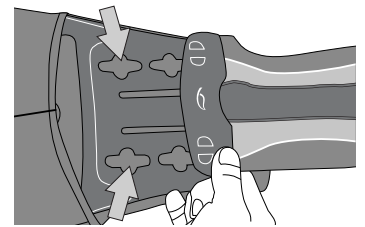

- C. Adjust length of tightening mechanism. For individual patient, it may be necessary to adjust length of closure string. Trim and adjust length of strings.

YES. AMOUNT CUT _____ NO

- D. If sizing yields extra plastic and if appropriate to individual's anatomy, trim extra plastic for superior customization to patient's individual anatomy.

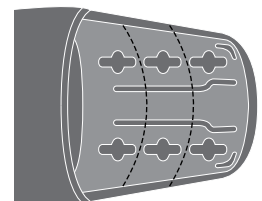

TIME SPENT: _____

VISTA® 637 LSO

Doctor: _____ Fitter: _____

Patient Name: _____ Date: _____



Patient #: _____ Additional Follow-Up Dates: _____

TOOLS NECESSARY: Scissors • Heat Gun • Tape Measure

STEP 4 - MODIFY RIGID PANELS

MODIFY ANTERIOR PANEL AND LATERAL PANELS AS NECESSARY

TIME SPENT: _____

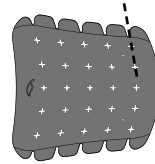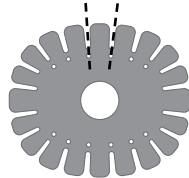


- Remove and trim to accommodate small and extra small anatomy.
- Remove and heat mold anterior panel as necessary.

STEP 5 - CUSTOMIZE BELT FIT


ANGLE ANTERIOR PANELS

Every patient has a unique individual anatomy. Determine angulation for proper fit. Circumferential contact at both upper and lower margins of brace is essential for proper brace performance and support.


- A. Bend anterior panel to conform to patient's anatomy.
- B. Angle anterior panels:

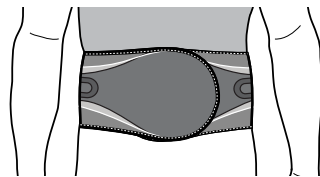

- Neutral Configuration for best support

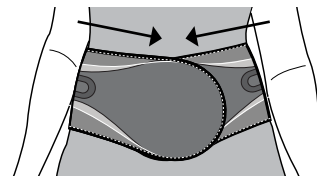

- Inferior Angulation Configuration for best support

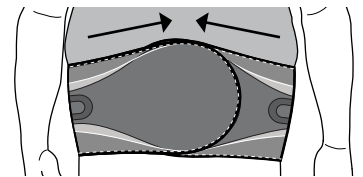

- Superior Angulation Configuration for best support

TIME SPENT: _____

STEP 6 - EDUCATION

EDUCATE PATIENTS

Proper education is needed for individual to maintain proper fit throughout total time of wear.

Items to educate patients on:

- Independent compression mechanics
- Don and doffing

- Proper angulation to ensure circumferential contact
- Proper placement of brace

- Proper cleaning
- Follow up appointments

TIME SPENT: _____

CLINICAL JUSTIFICATION FOR CUSTOMIZING BRACE

TOTAL TIME TO CUSTOMIZE BRACE: _____